

## HIGHLANDER® PROSTHETIC FOOT (FS3)

*Proven technology for dependable performance*

- Engineered for active amputees with long residual limbs
- Split toe design provides stability enhancing inversion and eversion
- Transfemoral amputees favor this foot due to its ideal blend of stability and action
- Exoskeletal configuration available upon request
- Heavy duty option is available for users weighing up to 500 lbs



HIGHLANDER® (FS3)

# HIGHLANDER® PROSTHETIC FOOT (FS3)

**Size\*: 22-30 cm**

\*Special Order Size 31 cm (7-10 day lead time)

**User Weight Rating:** 365 lbs

**Stiffness Category:** 1-9

**Connector:** Integral male pyramid


**Heel Height:** 3/8"

**Special Orders:**

Include request in the notes section on the order form.

- Exoskeletal Configuration (1 day lead time)
- HD Set-up for users weighing up to 500 lbs (7-10 day lead time)

**Warranty:** 36 months

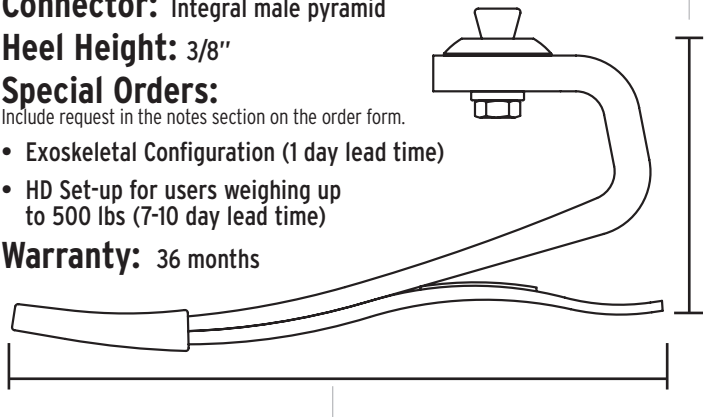

## Clearance & Weight:

Size	Clearance	Avg. Weights	
		Foot	Shell
22-25 cm	3.7 in	310g	127g
26-28 cm	4.2 in	371g	159g
29-31 cm	4.7 in	461g	219g

Note: A Spectra sock and pyramid connector are included with purchase.

## REQUIRED ACCESSORIES

### Cosmetic Foot Shell:

- High-profile (Default) FTC-2H
- Mid-profile (Optional) FTC-2M
- Choice of caucasian or ebony skin tone

## OPTIONAL ACCESSORIES

- **HD 34mm Component Kit, rated to 500lbs**  
(Pylon, 2 Female Tube Clamps, 4 Ear Connector) Part # KIT-00-202
- **Foot Shell Removal Tool-** Part #ACC-00-103
- **Additional Spectra Sock-** Part #ACC-00-101xx  
(choose 12", 14", or 18" length)
- **Additional Stiffening Bumpers-**
  - KIT-00-110 40 Durometer, Sizes 22-25cm
  - KIT-00-112 40 Durometer, Sizes 26-28cm
  - KIT-00-114 40 Durometer, Sizes 29-31cm