





The **Silhouette[®]** VS
Prosthetic Foot
offers curb appeal...

...with unquestionable **Shock Value**

Shock absorption and multi-axial ground compliance
soften otherwise unforgiving surfaces

Vertical shock is achieved through a unique foot and sole plate
system that facilitates displacement past the urethane pad
located at the heel.

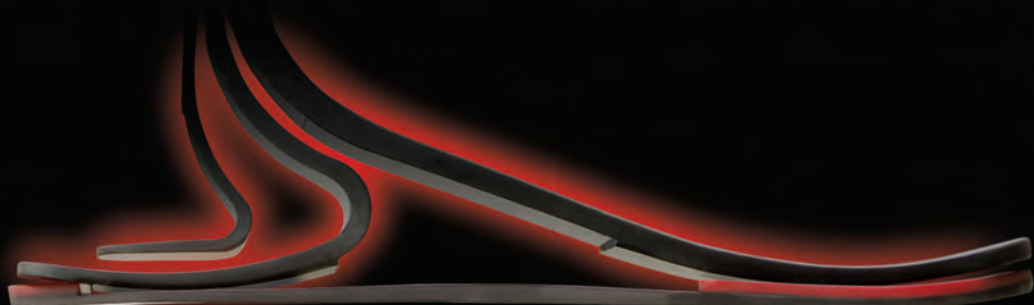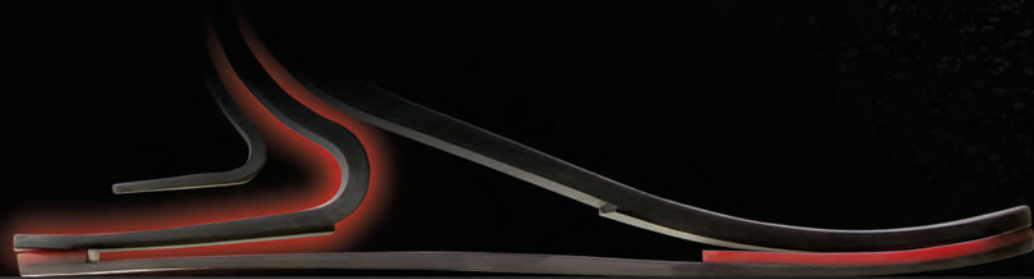

The Silhouette® VS

- **Vertical shock absorption softens impact at heel strike** to provide comfort and preserve residual limb health
- **Multi-axial function and +/-15 degrees of inversion/eversion** provide excellent ground compliance and stability
- The lightweight, slim-profile design facilitates an **easy to achieve cosmetic finish** capable of pleasing the most discerning critic
- The carbon fiber sole plate with tri-durometer urethane strips facilitates a **biomechanically correct ground reaction force**, moving lateral to medial through the foot with each step
- A **365 lb. weight rating** makes the Silhouette VS appropriate for a broad range of low to moderate impact K3 level ambulators
- The choice between **3/8" and 3/4" heel heights** provides the ability to wear higher heel dress shoes or work boots while maintaining proper prosthetic alignment



Inversion/eversion provides excellent ground compliance

An optimized sole plate facilitates shock absorption and a natural, medial roll-over



Silhouette® VS

Stiffness Category Selection Chart

Weight (lbs)	Category			Weight (kg)
	Low	Moderate	High	
100-115	1	1	2	44-52
116-130	1	2	3	53-59
131-150	2	3	4	60-68
151-170	3	4	5	69-77
171-195	4	5	6	78-88
196-220	5	6	7	89-100
221-255	6	7	8	101-116
256-285	7	8	9	117-130
286-325	8	9		131-147
326-365	9			148-166

Product Clearance and Weight

	Size	Clearance		Weight (Foot/Shell)	
		3/8 in Heel	3/4 in Heel		
Silhouette VS	22-25 cm	7.2 in	7.7 in	326g	127g
	26-28 cm	7.7 in	8.2 in	391g	159g
	29-31 cm	8.1 in	8.6 in	474g	219g

Part Number	Description	Heel Height
R14-00-xxxxx-xx	Silhouette VS Prosthetic Foot	3/8" heel
R18-00-xxxxx-xx	Silhouette VS Prosthetic Foot	3/4" heel

Specifications

Target user:

Low to moderate impact K3 level ambulators

Sizes: 22-30cm (31cm special order)

User Weight Rating: 365 lbs. (166kg)

Stiffness Categories: 1-9

Heel Height: 3/8" or 3/4"

Warranty: 36-month carbon fiber warranty

Connectors Kits (Choose One)

34mm female tube clamp (ACC-00-123 & ACC-00-124)

Male pyramid (ACC-00-196 & ACC-00-124)

Optional Accessories

Versa Clamp for pin systems and pylon lengthening (KIT-00-123V)

Knee Pyramid Adaptor: ACC-00-195

Suggested LCodes: L5987, L5986



30 Fairbanks • Ste. 114 • Irvine, CA 92618 • toll free: 888-818-6777 • ph: 949-672-0032
www.freedom-innovations.com

Also distributed by SPS, call 800-767-7776

© 2009 Freedom Innovations, LLC. All rights reserved. (11/09)