

# Fit, Finish and Roll.



Freedom DynAdapt™  
Carbon Fiber Foot

# Freedom DynAdapt™

## Carbon Fiber Foot

The Freedom DynAdapt™ foot is a slim profile, carbon fiber foot system with a slender, anatomic design for easy **fit** and **finish**. Its multi-axial function provides maximum comfort and the uninterrupted strands of carbon fiber in the full length heel provide patients with effortless **roll**over and a more natural gait.



### Fit, Finish...


Low profile, slender design results in an ankle that fits the anatomical design of the ankle, thus making the foot easy to cover cosmetically; and easy to fit & finish.

### ...and Roll

Full length attached heel with no bolts or bushings eliminates dead spots and weak structural areas, thus delivering a smooth rollover and comfortable walking with a more natural gait.

#### DynAdapt Force Deflection Data

Smooth line transitions of Force and Torque produced during a single step results in a more natural rollover performance.

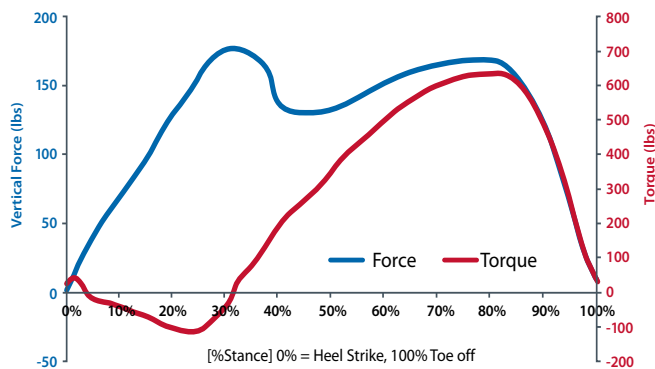

#### Blind Experiment (N=6)

Subjects who preferred the rollover performance of DynAdapt foot over competitive foot.



#### Maximum Comfort

Weighing in at just 380 g, the DynAdapt foot is more comfortable and requires less exertion



#### Increased Stability

Split toe/heel design provides excellent inversion and eversion for patient stability on uneven terrains

## Enhanced Stability


Split keel and heel design increase ground compliance and improve stability, thus allowing patients to walk with confidence on uneven surfaces or terrain.



### Reinforced for Strength


EnduraCore technology provides more durable and energy efficient performance.

- 1 Impact resistant, top surface film
- 2 High-compression fiber layer
- 3 High-modulus core technology
- 4 High-tension fiber layer
- 5 Impact resistant, bottom surface film

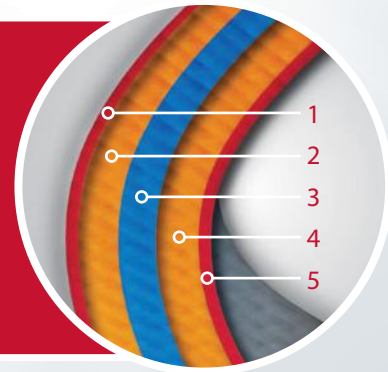

## Target Users

- All K3 level users who would benefit from improved stability and comfort
- Transtibial and transfemoral; Unilateral and bilateral amputees

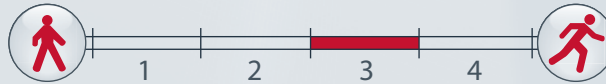

## Best All-Around Performance

- **For Users:** Combination of smooth rollover function, improved stability, and dynamic response makes DynAdapt the best all-around foot for all K3 users
- **For Prosthetists:** Ideal foot for prosthetists who are looking for a lightweight, durable, low profile graphite foot that is easy to cover cosmetically

### Smooth Rollover

A full length, unbolted sole plate ensures seamless stance phase roll-over for superior comfort

### Patient Trials (N=9)

Subjects who preferred the DynAdapt foot over their current foot.





Freedom Innovations is solely focused on developing world class lower limb solutions in close collaboration with amputees and prosthetists. These solutions encompass technological innovations, service, training, educational resources and consultative support to help ensure user satisfaction.



## Freedom DynAdapt™ Technical Specifications

**Sizes:** 22-31 cm

**Build Height:** 170 mm (Size 27)

**Product Weight:** 380 g (13.4 oz)  
(Size 27, Cat 5 graphite foot, without the footshell )

**User Weight Limit:** 166 kg (365 lbs)

**Stiffness Categories:** 9 (low, mod and high impact)

**Split Toe (inv/ev):** Yes

**Heel Height:** 10 mm (3/8 in)

**Foot Shell CAP:** Upon request

**Skin Tones:** 3

**Attachment Type:** Male pyramid

**Warranty:** 36 months (shell 6 months)

## Determine Your Impact Level

Low	Moderate	High
Walking	Jogging	Baseball
Gardening	Tennis	Football
Shopping	Hiking	Wakeboarding
Golfing	Volleyball	Snowboarding
Fishing	Skating	Soccer

## Stiffness Category Selection Chart

Weight (kg)	Impact Level			Weight (lbs)
	Low	Moderate	High	
44-52	1	1	2	100-115
53-59	1	2	3	116-130
60-68	2	3	4	131-150
69-77	3	4	5	151-170
78-88	4	5	6	171-195
89-100	5	6	7	196-220
101-116	6	7	8	221-255
117-130	7	8	9	256-285
131-147	8	9	Special Order	286-325
148-166	9	Special Order	Special Order	326-365

### Global Headquarters

3 Morgan  
Irvine, CA 92618  
U.S.A.  
toll-free 888.818.6777  
tel 949.672.0032  
fax 949.672.0084

[www.freedom-innovations.com](http://www.freedom-innovations.com)

### European Head Office

Jaargetijdenweg 4  
7532 SX Enschede  
The Netherlands  
toll-free 00800-2806 2806  
tel +31 (0)53-20 30 300  
fax +31 (0)53-20 30 305

[www.freedom-innovations.eu](http://www.freedom-innovations.eu)

