


Freedom Agilix™
Prosthetic Foot

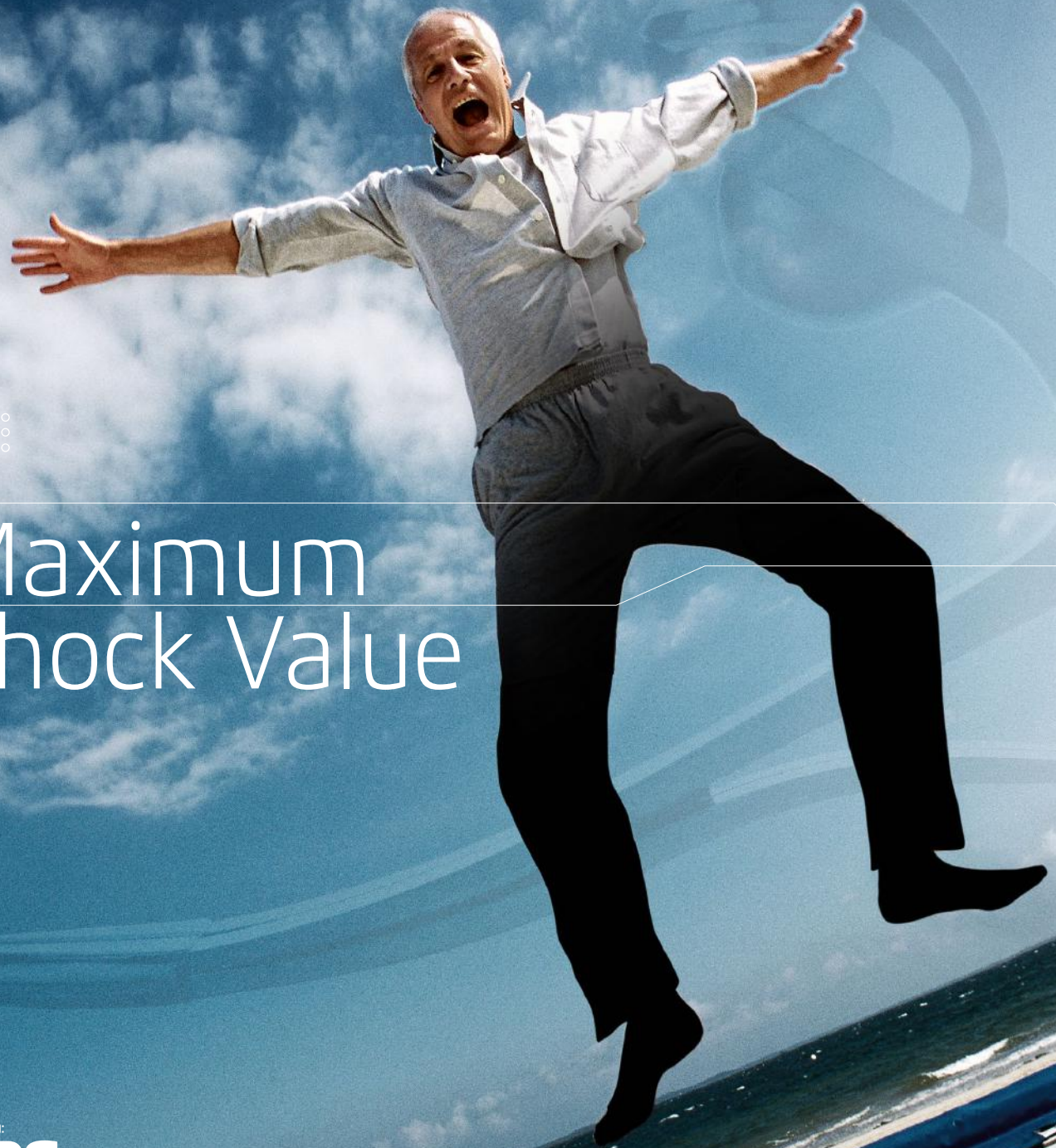

...

Maximum Shock Value

In Partnership with:

SPS

More than distribution.
Partnership.

 **FREEDOM**
INNOVATIONS



The Freedom Agilix™ is a multi-axial, shock absorbing flexible foot system designed to manage loading impacts, reduce socket shear forces and improve comfort while walking on nearly any terrain. The ultra-light-weight design offers K3 ambulators shockingly comfortable performance at a **value like no other.**

Freedom Agilix™ Prosthetic Foot

Multi-Axial Shock Absorbing
Foot System



Target Users

K3 level users who would benefit from increased shock absorption, stability and comfort. Rated for users weighing up to 365 lbs.

- Users who prefer a softer, easy to roll-over, flexible foot system
- Users who are sensitive to socket pressures and prone to limb discomfort
- Users who walk for extended periods for their vocation or exercise



Not recommended: High impact or K4 level users; athletes without a secondary prosthesis



EVERYDAY

Freedom Agilix is part of the "Everyday Prosthetic Foot Series"

How it Works

Shock Absorption



Shock absorption and multi-axial function in a sleek, integrated design make the Agilix lightweight and reliable

Vertical shock absorption and carbon fiber flexibility optimized for low to moderate impact K3 ambulators; reduces socket shear forces and improves overall comfort

Multi-axial function

increases ground compliance and stability; Streamlined, split-toe design results in improved long-term durability, saving time and money

Full length dynamic heel attached at the toe ensures a seamless roll-over, for superior comfort and a more symmetrical gait

Lower profile and lighter weight than the competitive alternative, allowing users to do more and feel less tired at the end of the day



Focused Solutions

Freedom Innovations is solely focused on developing world-class lower-limb solutions in close collaboration with amputees and prosthetists. These solutions encompass technological innovations, service, training, educational resources and consultative support to help ensure user satisfaction.

Freedom Agilix™ Technical Specifications

Sizes: 22-31 cm

(Sandal-toe available for sizes 22-28cm)

Build Height: 161 mm (Size 27)

Product Weight: 650 g/ 23 oz
(27cm, cat 4, incl. shell)

User Weight Limit: 166 kg (365 lbs)

Stiffness Categories: 9 (low and mod impact)

Split Toe (inv/ev): Yes

Heel Height: 10 mm ($\frac{3}{8}$ in)

Foot Shell CAP: Upon request

Skin Tones: 3

Attachment Type: Male pyramid

Warranty: 36 months (shell 6 months)

Stiffness Category Selection Chart

Weight (lbs)	Impact Level		Weight (kg)
	Low	Moderate	
100-115	1	1	44-52
116-130	1	2	53-59
131-150	2	3	60-68
151-170	3	4	69-77
171-195	4	5	78-88
196-220	5	6	89-100
221-255	6	7	101-116
256-285	7	8	117-130
286-325	8	9	131-147
326-365	9		148-166



Determine Your Impact Level

Low	Moderate
Walking	Light jogging
Gardening	Tennis
Bowling	Volleyball
Shopping	Hiking
Golfing	Dancing
Biking	Skating
Fishing	
Camping	

Global Headquarters

3 Morgan
Irvine, CA 92618
toll free: 888-818-6777 phone: 949-672-0032
fax: 949-672-0084
www.freedom-innovations.com

European Head Office

Jaargetijdenweg 4
7532 SX Enschede
The Netherlands
toll free: 00800 - 2806 2806 phone: 31 (0) 53 20 30 300
fax: 31 (0) 53 20 30 305
www.freedom-innovations.eu

Manufacturing & Returns

Freedom Innovations
425 East 400 North
Gunnison, UT 84634

In Partnership with:
SPS
More than distribution.
Partnership.

