

# Silhouette<sup>®</sup> VS & Silhouette<sup>®</sup> LP-VS

Shock absorption and multi-axial ground compliance  
soften otherwise unforgiving surfaces

- Vertical shock absorption softens impact at heel strike to provide comfort and preserve residual limb health
- Multi-axial function and +/- 15 degrees of inversion/eversion provide excellent ground compliance and stability
- The lightweight, slim-profile design facilitates an easy to achieve cosmetic finish capable of pleasing the most discerning critic
- The carbon fiber sole plate with tri-durometer urethane strips facilitates a biomechanically correct ground reaction force, moving lateral to medial through the foot with each step
- Available in 10 mm (3/8in) heel height
- A 166 kg (365 lbs) weight rating makes the Silhouette VS appropriate for a broad range of low to moderate impact K3 level ambulators
- 36-month warranty



Silhouette VS (R14)



Silhouette LP-VS (R15)



Vertical Shock



Multi-Axial



STYLE

**Product ID:** SilhouetteVS (R14);

SilhouetteLP-VS (R15)

**Sizes:** 22-31 cm


**User Weight Rating:** 166 kg (365 lbs)

**Heel Height:** 10mm (3/8")

**Connector:** Male or female pyramid, or Versa Clamp

**Warranty:** Foot module 36 months,

Cosmetic foot shell 6 months

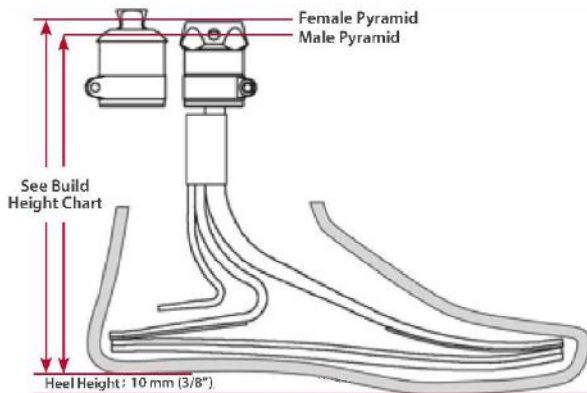

Note: Build height measurements represent a fully assembled unit.

Note: Build height measurements represent a fully assembled unit.

## Accessories

- Female pyramid, 34 mm: (ACC-00-129)
- Male pyramid, 34 mm: (ACC-00-196)
- Tube Clamp Insert: (ACC-00-124)  
(custom cut, foot module serial number required)
- Versa Clamp Adapter Kit: (KIT-00-123V) includes extended female pyramid and spacer ring
  - Extended 34 mm Female TubeClamp: (ACC-EX-123)
  - Spacer Ring: (KIT-00-123S)
- Knee Pyramid Adapter, Male: (ACC-00-195) for knees receiving a 30 mm pylons
- Replacement Sole Plate, Tri Durometer: (RS4-SP-0xAxx-xx)
- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx) 12", 14", 16"
- Stiffening Bumpers:
  - 22-25 cm (KIT-00-12101-00)
  - 26-28 cm (KIT-00-12202-00)
  - 29-31 cm (KIT-00-12203-00)

## Build Height & Weight

|  | Size (cm) | Height (in / mm) | Weight (oz / g) |
|--|-----------|------------------|-----------------|
| <b>Silhouette VS w/female adapter</b>  | 22-25     | 7 / 178          | 19 / 525        |
|  | 26-28     | 7½ / 191         | 23 / 645        |
|  | 29-31     | 7¾ / 201         | 28 / 785        |
| <b>Silhouette LP-VS w/male adapter</b> | 22-25     | 6¾ / 168         | 15 / 430        |
|  | 26-28     | 7⅛ / 181         | 19 / 550        |
|  | 29-31     | 7½ / 191         | 24 / 685        |

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

## Stiffness Category Selection

| Weight (lbs) | Impact Level |          |      | Weight (kg) |
|--------------|--------------|----------|------|-------------|
|              | Low          | Moderate | High |             |
| 100-115      | 1            | 1        | 2    | 44-52       |
| 116-130      | 1            | 2        | 3    | 53-59       |
| 131-150      | 2            | 3        | 4    | 60-68       |
| 151-170      | 3            | 4        | 5    | 69-77       |
| 171-195      | 4            | 5        | 6    | 78-88       |
| 196-220      | 5            | 6        | 7    | 89-100      |
| 221-255      | 6            | 7        | 8    | 101-116     |
| 256-285      | 7            | 8        | 9    | 117-130     |
| 286-325      | 8            | 9        |      | 131-147     |
| 326-365      | 9            |          |      | 148-166     |

## Foot Module Part Number

|   |   |                        |                      |   |   |                      |   |                      |
|---|---|------------------------|----------------------|---|---|----------------------|---|----------------------|
| <input type="text"/>                            | -0  | <input type="text"/>   | -0                   | <input type="text"/>  | A | <input type="text"/> | - | <input type="text"/> |
| <b>Product ID</b><br>R14 (10 mm)<br>R15 (19 mm) | <b>Adapter</b><br>1-female<br>2-male<br>3-versa | <b>Category</b><br>1-9 | <b>Size</b><br>22-31 | <b>Toe Configuration</b><br>RR-regular toe, right<br>RL-regular toe, left |   |                      |   |                      |

## Foot Shell Part Number

|                          |                                   |                             |                      |   |                                  |                      |                      |                      |     |                      |
|--------------------------|-----------------------------------|-----------------------------|----------------------|---|----------------------------------|----------------------|----------------------|----------------------|-----|----------------------|
| <input type="text"/>     | FTC-                              | <input type="text"/>        | <input type="text"/> | -   | <input type="text"/>             | <input type="text"/> | <input type="text"/> | <input type="text"/> | 4-R | <input type="text"/> |
| <b>Generation</b><br>2-3 | <b>Profile</b><br>M-mid<br>H-high | <b>Cap</b><br>0-yes<br>1-no | <b>Size</b><br>22-31 | <b>Skin Tone</b><br>L-light<br>M-medium<br>D-dark | <b>Side</b><br>L-left<br>R-right |                      |                      |                      |     |                      |

See the Cosmetic Foot Shell page for more detailed ordering instructions.

