


Silhouette® & Silhouette LP®

High science meets high style with simplified fitting

- Strong, lightweight carbon-fiber construction delivers an extended range of motion allowing user to walk farther, faster, using less energy
- Slim profile facilitates an easy and refined cosmetic finish
- Advanced dual heel system provides a dynamic response for normal ambulation and an override spring for security during higher activities
- Ideal for low to moderate impact K3 amputees
- Optional low profile version with an integral male pyramid simplifies the fitting process
- Versa Clamp provides an easy fix for pylons cut too short or delivers the space needed for a locking pin (sold separately)
- Available in 10 mm (3/8") heel height
- Rated up to 166 kg (365 lbs)
- 36-month warranty

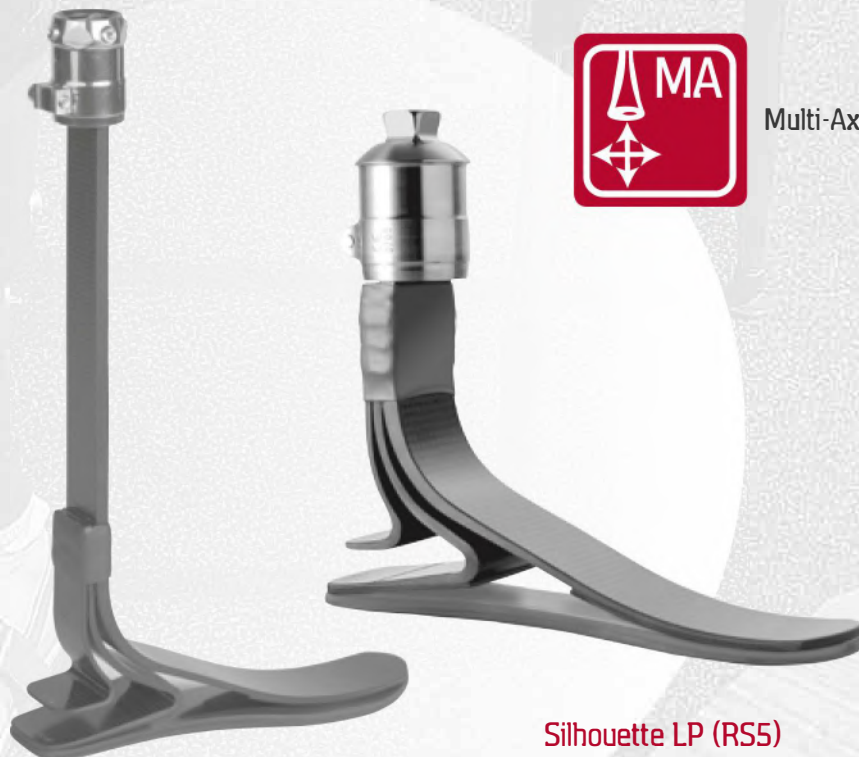

Silhouette (RS4)

Silhouette LP (RS5)



STYLE

Product ID: Silhouette (RS4); Silhouette LP (RS5)


Sizes: 22-31 cm

User Weight Rating: 166 kg (365 lbs)

Heel Height: 10 mm (3/8")

Connector: Male or female pyramid, or Versa Clamp

Warranty: Foot module 36 months,
Cosmetic foot shell 6 months

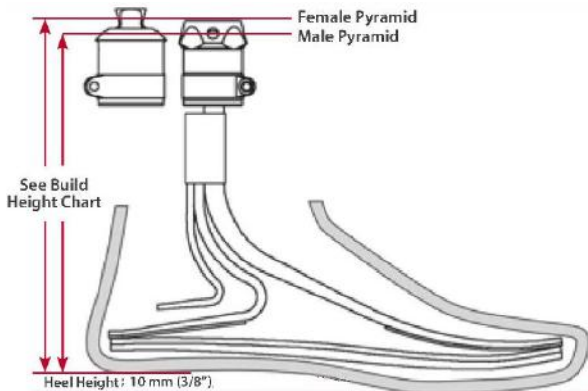

Note: Build height measurements represent a fully assembled unit.

Note: Build height measurements represent a fully assembled unit.

Accessories

- Female pyramid, 34 mm: (ACC-00-129)
- Male pyramid, 34 mm: (ACC-00-196)
- Tube Clamp Insert: (ACC-00-124)
(custom cut, foot module serial number required)
- Versa Clamp Adapter Kit: (KIT-00-123V) includes extended female pyramid and spacer ring
 - Extended 34 mm Female TubeClamp: (ACC-EX-123)
 - Spacer Ring: (KIT-00-123S)
- Knee Pyramid Adapter, Male: (ACC-00-195) for knees receiving a 30 mm pylons
- Replacement Sole Plate, Tri Durometer: (RS4-SP-0xAxx-xx)
- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx) 12", 14", 16"
- Stiffening Bumpers:
 - 22-25 cm (KIT-00-12101-00)
 - 26-28 cm (KIT-00-12202-00)
 - 29-31 cm (KIT-00-12203-00)

Build Height & Weight

	Size (cm)	Height (in / mm)	Weight (oz / g)
Silhouette w/female adapter	22-25	7 / 178	19 / 525
	26-28	7½ / 191	23 / 645
	29-31	7¾ / 201	28 / 785
Silhouette LP w/male adapter	22-25	6½ / 168	15 / 430
	26-28	7½ / 181	19 / 550
	29-31	7½ / 191	24 / 685

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

Stiffness Category Selection

	Impact Level			
	Low	Moderate	High	
100-115	1	1	2	44-52
116-130	1	2	3	53-59
131-150	2	3	4	60-68
151-170	3	4	5	69-77
171-195	4	5	6	78-88
196-220	5	6	7	89-100
221-255	6	7	8	101-116
256-285	7	8	9	117-130
286-325	8	9		131-147
326-365	9			148-166

Foot Module Part Number

Product ID	Adapter	Category	Size	Toe Configuration
RS4 (10 mm) RSS (LP 10 mm)	1 female 2-male 3-versa	1-9	22-31	RR-regular toe, right RL-regular toe, left

Foot Shell Part Number

Generation	Profile	Cap	Size	Skin Tone	Side
2-3	M-mid H-high	0 yes 1-no	22-31	L-light M-medium D-dark	L-left R-right

See the Cosmetic Foot Shell page for more detailed ordering instructions.

