

Thrive[®]

The world's first load-activated prosthetic foot

- Responds automatically when lifting or carrying heavy objects
- A unique dual-keel design utilizes a full-length primary keel and a secondary, load-activated keel to ensure consistent response
- Full-length carbon fiber heel plate technology facilitates more natural ambulation with and without variable loads
- Responds to variable loads of up to 30% of the user's body weight
- Split toe construction with +/-15 degrees of inversion and eversion creates dynamic control on irregular surfaces
- Carbon fiber and sole plate design provides 11 mm of vertical compression for increased control and comfort
- Rated to 166 kg (365 lbs)
- 36-month warranty



Thrive (FS5)



Variable Stiffness



Split Toe



Sandal Toe



ACTIVE

Product ID: FS5

Sizes: 22-31 cm


Sandal-toe available for sizes 22-28 cm

User Weight Rating: 166 kg (365 lbs)

Heel Height: 10 mm (3/8")

Connector: Male pyramid

Warranty: Foot module 36 months,
Cosmetic foot shell 6 months

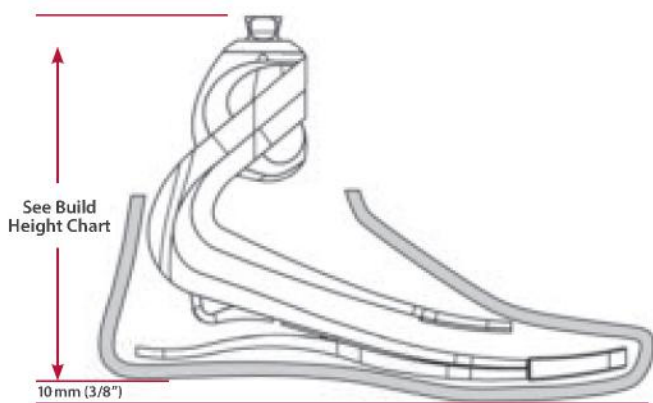

Note: Build height measurements represent a fully assembled unit.

Accessories

- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Socks:
 - 22-25 cm: (KIT-00-FS525-0C)
 - 26-28 cm: (KIT-00-FS528-0D)
 - 29-31 cm: (KIT-00-FS531-0E)
- Stiffening Bumpers:
 - 22-25 cm: (KIT-00-110)
 - 26-28 cm: (KIT-00-112)
 - 29-31 cm: (KIT-00-114)

Build Height & Weight

Size (cm)	Height (in / mm)	Weight (oz / g)
22-25	5¾ / 145	22 / 630
26-28	6½ / 157	27 / 760
29-31	7 / 175	31 / 870

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

Stiffness Category Selection

Weight (lbs)	Impact Level			Weight (kg)
	Low	Moderate	High	
100-115	1	1	2	44-52
116-130	1	2	3	53-59
131-150	2	3	4	60-68
151-170	3	4	5	69-77
171-195	4	5	6	78-88
196-220	5	6	7	89-100
221-255	6	7	8	101-116
256-285	7	8	9	117-130
286-325	8	9		131-147
326-365	9			148-166

Foot Module Part Number

FS5-00-0 **A** -

Category: 1-9 Size: 22-31 Toe Configuration: RU-regular toe, universal; SR-sandal toe, right; SL-sandal toe, left

Foot Shell Part Number

FTC- - **4-**

Generation: 2-3 Profile: M-mid, H-high Cap: 0-yes, 1-no Size: 22-31 Skin Tone: L-light, M-medium, D-dark Toe: R-regular, S-sandal Side: L-left, R-right

See the Cosmetic Foot Shell page for more detailed ordering instructions.

FEET

