

K3/K4

FREEDOM

MAVERICK

TREME



Patient Selection Guide

Profile of **Ideal User**

The Freedom Innovations' Maverick Xtreme foot was designed for active K3 and K4 level ambulators that require a product with increased durability without compromising the desired energy storage/return and smooth roll-over during stance phase of gait. The increased durability is a result of the fiberglass material used to manufacture the foot and a result of the innovative heel/keel design. The material combined with the design allow for a natural, flexible feel during activities of daily living (ADL's) and various strenuous exercises. This new lighter weight fiberglass option is waterproof and is available in a sandal toe option for greater lifestyle versatility.

Features

- Coding: L5987 / L5986
- High impact, extreme user
- Regular and Sandal Toe Option
- Waterproof
- Available in sizes 23cm - 30cm
- Light weight fiber glass option
- Increased durability, flexibility and natural feel
- Proudly Made in the USA

Specifications

- User Weight Ratings - 365 lbs.
- Build Height - 6.8" / 175mm
- Average Product Weight*: 1.63lbs / 743g (with footshell)
- Warranty Foot Module 36 months, cosmetic footshell 6 months



*Based on size 26cm, category 4

FREEDOM
INNOVATIONS

Ordering Guide

1 Verify User's Measurements

- Patient weight
- Foot Size (23-30cm)
- Right or Left
- Sandal Toe*/ Standard Toe


*Available Sandal Toe Sizes: 23cm-28cm

2 Select Stiffness Category



Weight (lbs)	Impact Level			Weight (kg)
	Low	Moderate	High	
100-115	1	1	2	44-52
116-130	1	2	3	53-59
131-150	2	3	4	60-68
151-170	3	4	5	69-77
171-195	4	5	6	78-88
196-220	5	6	7	89-100
221-255	6	7	8	101-116
256-285	7	8	9	117-130
286-325	8	9		131-147
326-365	9			148-166

3 Complete Foot Module Part Number

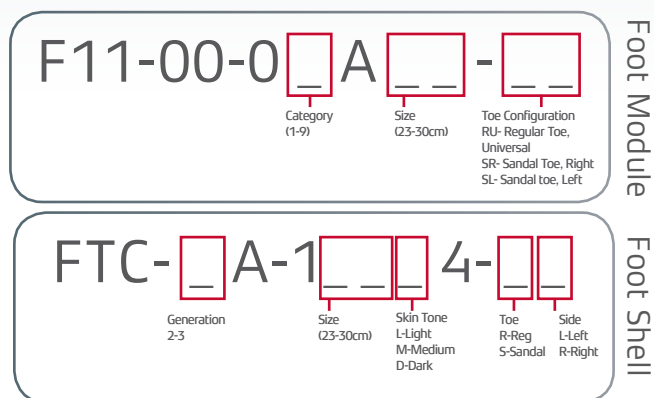

Global Headquarters
 3 Morgan
 Irvine, CA 92618
 toll free: 888-818-6777
 phone: 949-672-0032
 fax: 949-672-0084
www.freedom-innovations.com

European Head Office
 Jaargetijdenweg 4
 7532 SX Enschede
 The Netherlands
 toll free: 00800 - 2806 2806
 phone: 31 (0) 53 20 30 300
 fax: 31 (0) 53 20 30 305
www.freedom-innovations.eu

USA Manufacturing & Returns
 425 East 400 North
 Gunnison, UT 84634

Connect with us



#freedommaverick