

Renegade® A·T & Renegade® LP-A·T

All·Terrain Performance AMPLified

- The most dynamic, versatile and durable Renegade prosthetic foot
- EnduraCore® Technology delivers 3 times the fatigue life of standard laminates
- Unique design returns up to 35 % more energy than competitive alternatives
- Split-toe design delivers +/- 15 degrees of inversion/eversion for enhanced stability
- Z-Shock Technology™ absorbs shear forces providing greater user comfort
- Lightweight construction – only 515 grams
- Rated to 166 kg
- 36-month warranty



Renegade A·T (R11)

Renegade LP-A·T (R16)



Vertical Shock




Split Toe



EnduraCore® Technology

Renegade® A·T & Renegade® LP-A·T

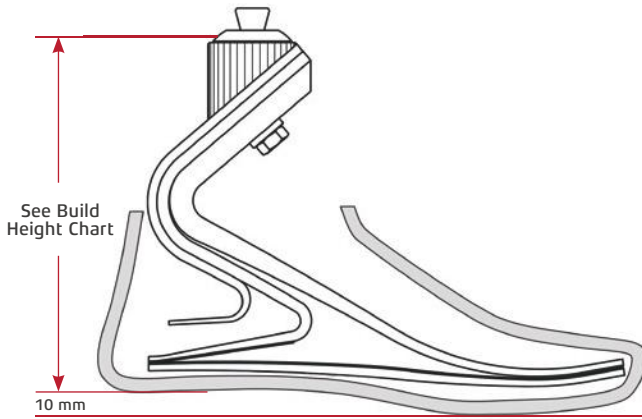

Note: Build height measurements represent a fully assembled unit.

Product ID: Renegade A·T (R11); Renegade LP-A·T (R16)

Sizes: 22-31 cm

User Weight Rating: 166 kg

Heel Height: 10 mm

Connector: Male pyramid

Warranty: Foot module 36 months,
Cosmetic foot shell 6 months

Accessories

- Pre-Cut Shin Foam: (R11 only; two piece clamshell design)
 - 22-25 cm (AFP-00-RS101)
 - 26-28 cm (AFP-00-RS102)
 - 29-31 cm (AFP-00-RS103)
- Sole Plate: (size and category specific)
 - Universal toe: (RS1-SP)
- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx)
 - 30 cm, 35 cm, 40 cm
- Stiffening Bumpers:
 - 22-25 cm (KIT-00-116)
 - 26-28 cm (KIT-00-117)
 - 29-31 cm (KIT-00-118)

Note: Spectra socks and heel bumpers are included in the delivery.

Build Height & Weight

	Size (cm)	Height (mm)	Weight (g)
Renegade A·T	22-25	162	625
	26-28	184	765
	29-31	200	945
Renegade LP-A·T	22-25	127	625
	26-28	140	760
	29-31	149	930

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

Stiffness Category Selection

Impact Level			Weight (kg)
Low	Moderate	High	
1	1	2	44-52
1	2	3	53-59
2	3	4	60-68
3	4	5	69-77
4	5	6	78-88
5	6	7	89-100
6	7	8	101-116
7	8	9	117-130
8	9		131-147
9			148-166

Complete Prosthetic Foot (incl. foot shell) Part Number

-00-0A-M-3

Product ID: R11 (R16)
 Category: 1-9
 Size: 22-31
 Side: L-left (R-right)
 Skin Tone: L-light (M-medium, D-dark)

Foot Module Part Number

-00-0A-RU

Product ID: R11 (Renegade A·T), R16 (Renegade LP-A·T)
 Category: 1-9
 Size: 22-31

Foot Shell Part Number

FTC--4-R

Generation: 2-3
 Profile: M-mid (H-high)
 Cap: 0-yes (1-no)
 Size: 22-31
 Skin Tone: L-light (M-medium, D-dark)
 Side: L-left (R-right)

See the Cosmetic Foot Shell page for more detailed ordering instructions.