

SPECIALTY

# Runway® & Runway® HX

## Advanced technology delivers unmatched versatility

- As the only user-adjustable heel height foot designed for active K3 amputees, the device is customized to user weight and activity level to provide superior comfort and energy return
- Adjust heel height from zero to 5 cm; option of push button or 4 mm hex key adjustable design
- Push button design offers 10 height settings
- Runway HX offers infinite adjustments on the ankle slide
- Anatomical 'gliding' ankle design maintains appropriate prosthetic alignment, regardless of heel height
- Sandal toe foot shell offers improved user acceptance
- Rated to 116 kg
- 36-month carbon fiber warranty
- 18-month ankle mechanism warranty




Runway HX (RS2-K1)

Runway (RS2)




Sandal Toe

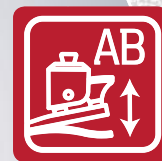

Height Adjustable

EnduraCore®  
Technology

# Runway® & Runway® HX

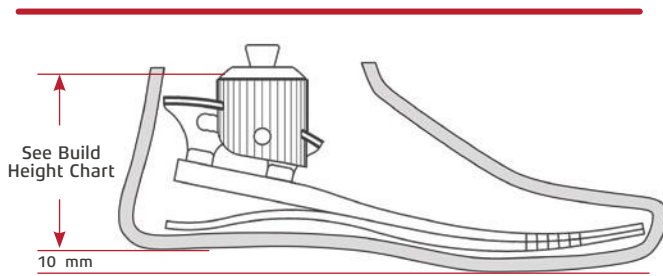

Note: Build height measurements represent a fully assembled unit.

**Product ID:** Runway (RS2); Runway HX (RS2-K1)

**Sizes:** 22-28 cm\*

**User Weight Rating:** 116 kg

**Heel Height:** 0-50 mm

**Connector:** Male pyramid

**Warranty:** Foot module 6 months, Ankle 18 months, Cosmetic foot shell 6 months

\*When ordering the Sandal Toe option, please contact Customer Service for lead time.

## Accessories

- Replacement Ankle: (ACC-00-10500-Rx)  
Left or Right
- Shim Adjustment Wrench: (ACC-00-104)
- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx)  
30 cm, 35 cm, 40 cm
- Stiffening Bumpers:
  - 22-25 cm (KIT-00-11313)
  - 26-28 cm (KIT-00-11314)

Note: Spectra socks and heel bumpers are included in the delivery.

## Build Height & Weight

Size (cm)	Height (mm)	Weight (g)
22-25	98	690
26-28	100	792

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

## Stiffness Category Selection

Impact Level			Weight (kg)
Low	Moderate	High	
1	1	2	44-52
1	2	3	53-59
2	3	4	60-68
3	4	5	69-77
4	5	6	78-88
5	6		89-100
6			101-116

## Complete Prosthetic Foot (incl. foot shell) Part Number

**RS2--1A---**

Product ID: 00 (Runway (button)), K1 (Runway HX (hex key))  
 Category: 1-6  
 Size: 22-28  
 Version: M-mid profile, S-sandal toe  
 Side: L-left, R-right  
 Generation: 3-regular toe, 2-sandal toe  
 Skin Tone: L-light, M-medium, D-dark

**Foot Module Part Number**

**RS2--0A-**

Product ID: 00 (Runway), K1 (Runway HX)  
 Category: 1-6  
 Size: 22-28  
 Toe Configuration: RR-regular toe, right; RL-regular toe, left; SR-sandal toe, right; SL-sandal toe, left

**Foot Shell Part Number**

**FTC-----4--**

Generation: 2-3  
 Profile: M-mid, H-high  
 Cap: 0=yes, 1=no  
 Size: 22-28  
 Skin Tone: L-light, M-medium, D-dark  
 Toe: S-sandal  
 Side: R-right

See the Cosmetic Foot Shell page for more detailed ordering instructions.