

SPECIALTY

# Silhouette® & Silhouette® LP

## High science meets high style with simplified fitting


- Strong, lightweight carbon fiber construction delivers an extended range of motion allowing user to walk farther, faster, using less energy
- Slim profile facilitates an easy and refined cosmetic finish
- Advanced dual heel system provides a dynamic response for normal ambulation and an override spring for security during higher activities
- Ideal for low to moderate impact K2 and K3 amputees
- Optional low profile version with an integral male pyramid simplifies the fitting process
- Versa Clamp provides an easy fix for pylons cut too short or delivers the space needed for a locking pin (sold separately)
- The choice between 10 mm and 19 mm heel heights
- Rated up to 166 kg
- 36-month warranty



# Silhouette® & Silhouette® LP



SPECIALTY

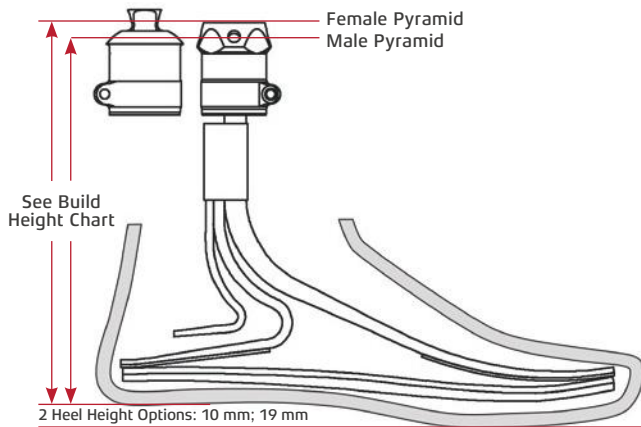

Note: Build height measurements represent a fully assembled unit.

**Product ID:** Silhouette (RS4/RS8);  
Silhouette LP (RS5/RS9)

**Sizes:** 22-31 cm

**User Weight Rating:** 166 kg

**Heel Height:** 10 mm or 19 mm

**Connector:** Male or female pyramid, or Versa Clamp

**Warranty:** Foot module 36 months,  
Cosmetic foot shell 6 months

## Accessories

- Female pyramid, 34 mm: (ACC-00-129)
- Male pyramid, 34 mm: (ACC-00-196)
- Tube Clamp Insert: (ACC-00-124)  
(custom cut, foot module serial number required)
- Versa Clamp Adapter Kit: (KIT-00-123V) includes extended female pyramid and spacer ring
  - Extended 34 mm Female TubeClamp: (ACC-EX-123)
  - Spacer Ring: (KIT-00-123S)
- Knee Pyramid Adapter, Male: (ACC-00-195) for knees receiving a 30 mm pylons
- Replacement Sole Plate, Tri Durometer: (RS4-SP-0xAxx-xx)
- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx)  
30 cm, 35 cm, 40 cm
- Stiffening Bumpers:
  - 22-25 cm (KIT-00-12101)
  - 26-28 cm (KIT-00-12102)
  - 29-31 cm (KIT-00-12103)

Note: Spectra socks and heel bumpers are included in the delivery.

## Build Height & Weight

	Size (cm)	Height (mm)	Weight (g)
Silhouette w/female adapter	22-25	178	525
	26-28	191	645
	29-31	201	785
Silhouette LP w/male adapter	22-25	168	430
	26-28	181	550
	29-31	191	685

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

## Stiffness Category Selection

Impact Level			Weight (kg)
Low	Moderate	High	
1	1	2	44-52
1	2	3	53-59
2	3	4	60-68
3	4	5	69-77
4	5	6	78-88
5	6	7	89-100
6	7	8	101-116
7	8	9	117-130
8	9		131-147
9			148-166

## Complete Prosthetic Foot (incl. foot shell) Part Number

---A-M-

<b>Product ID</b>	<b>Adapter</b>	<b>Category</b>	<b>Size</b>	<b>Side</b>	<b>Skin Tone</b>
RS4 (10 mm)	1-female	1-9	22-31	L-left	L-light
RS5 (LP 10 mm)	2-male			R-right	M-medium
RS8 (19 mm)	3-versa				D-dark
RS9 (LP 19 mm)					

### Foot Module Part Number

---A-

<b>Product ID</b>	<b>Adapter</b>	<b>Category</b>	<b>Size</b>	<b>Toe Configuration</b>
RS4 (10 mm)	1-female	1-9	22-31	RR-regular toe, right
RS5 (LP 10 mm)	2-male			RL-regular toe, left
RS8 (19 mm)	3-versa			
RS9 (LP 19 mm)				

### Foot Shell Part Number

**FTC-****4-R**

<b>Generation</b>	<b>Profile</b>	<b>Cap</b>	<b>Size</b>	<b>Skin Tone</b>	<b>Side</b>
2-3	M-mid H-high	0-yes 1-no	22-31	L-light M-medium D-dark	L-left R-right

See the Cosmetic Foot Shell page for more detailed ordering instructions.