



1

2

3

4



SPORT

# Slalom™ Ski

## Just step in and go!

- The first prosthetic foot designed to be quickly and easily inserted directly into a standard ski binding
- Eliminates the need for a ski boot
- Provides necessary flex and shock absorption for fluid motion
- Better control and range of motion increases endurance, allowing more time on the slopes
- Rated to 114 kg
- 12-month warranty, custom product




Slalom Ski (SKI)

# Slalom™ Ski



SPORT

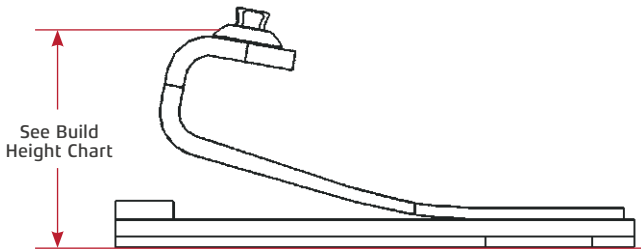

See Build Height Chart

Note: Build height measurements represent a fully assembled unit.

**Product ID:** SKI

**Sizes:** Custom, 25 cm minimum

**User Weight Rating:** 114 kg

**Connector:** Male pyramid

**Lead Time:** 10-14 days

**Warranty:** 12 months

**Return Policy:** No trial period

## How to Measure for Size

### Unilateral user

- Preferred measurement method: Measure the length of the user's sound limb ski boot in cm, from heel to toe
- Alternative method: Measure the length of the ski binding in cm

### Bilateral user

- Measure the length of the ski binding in cm

### Note

- The Slalom Ski Foot is made to “DIN specifications” in order to pair with standard ski bindings. If the user is using something other than standard ski bindings, please contact Customer Service to discuss special needs.

## Build Height & Weight

Size (cm)	Height (mm)	Weight (g)
Custom	140	1055

Measurements are approximate and will vary due to category selection and manufacturing variability.

## Stiffness Category Selection

Category	Weight (kg)
11	44-72
12	73-90
13	91-114

## Foot Module Part Number

SKI-00-A-RU

Category 11-13      Size cm (25 cm minimum)