



1

2

3

4



SPECIALTY

Thrive[®]

The world's first load-activated prosthetic foot

- Responds automatically when lifting or carrying heavy objects
- A unique dual-keel design utilizes a full-length primary keel and a secondary, load-activated keel to ensure consistent response
- Full-length carbon fiber heel plate technology facilitates more natural ambulation with and without variable loads
- Responds to variable loads of up to 30 % of the user's body weight
- Split toe construction with +/-15 degrees of inversion and eversion creates dynamic control on irregular surfaces
- Carbon fiber and sole plate design provides 11 mm of vertical compression for increased control and comfort
- Rated to 166 kg
- 36-month warranty



Thrive (FS5)




Split Toe



Sandal Toe

Thrive[®]

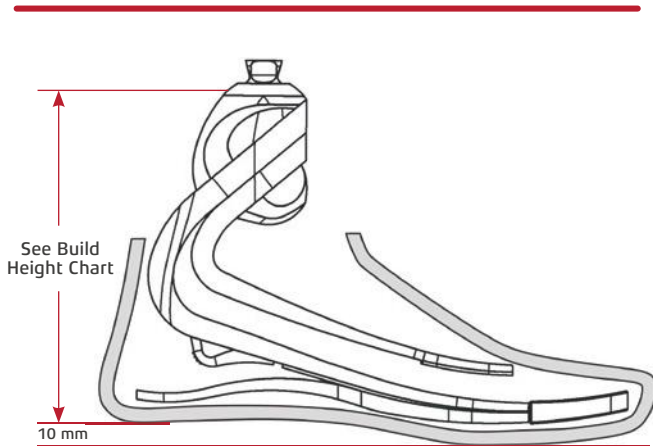

Note: Build height measurements represent a fully assembled unit.

Product ID: FS5

Sizes: 22-31 cm

Sandal toe available for sizes 22-28 cm*

User Weight Rating: 166 kg

Heel Height: 10 mm

Connector: Male pyramid

Warranty: Foot module 36 months,
Cosmetic foot shell 6 months

*When ordering the Sandal Toe option, please contact Customer Service for lead time.

Accessories

- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx)
30 cm, 35 cm, 40 cm
- Stiffening Bumpers:
 - 22-25cm (KIT-00-14101-00)
 - 26-28cm (KIT-00-14102-00)
 - 29-31cm (KIT-00-14103-00)

Note: Spectra socks and heel bumpers are included in the delivery.

Build Height & Weight

Size (cm)	Height (mm)	Weight (g)
22-25	145	630
26-28	157	760
29-31	175	870

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

Stiffness Category Selection

Impact Level			Weight (kg)
Low	Moderate	High	
1	1	2	44-52
1	2	3	53-59
2	3	4	60-68
3	4	5	69-77
4	5	6	78-88
5	6	7	89-100
6	7	8	101-116
7	8	9	117-130
8	9		131-147
9			148-166

Complete Prosthetic Foot (incl. foot shell) Part Number

FS5-00-A--

Toe	Category	Size	Version	Side	Generation	Skin Tone
0-regular 1-sandal	1-9	22-31	M-mid profile S-sandal toe	L-left R-right	3-regular toe 2-sandal toe	L-light M-medium D-dark

Foot Module Part Number

FS5-00-0A

Category	Size	Toe Configuration
1-9	22-31	RU-regular toe, universal SR-sandal toe, right SL-sandal toe, left

Foot Shell Part Number

FTC--4-

Generation	Profile	Cap	Size	Skin Tone	Toe	Side
2-3	M-mid H-high	0-yes 1-no	22-31	L-light M-medium D-dark	R-regular S-sandal	L-left R-right

See the Cosmetic Foot Shell page for more detailed ordering instructions.